

respiralogics™

Adding bubbles
to CPAP is proving
beneficial in treating
infants. Now you
have the perfect
complement of
products...

baby™ and

Babi.Plus®



FEATURES AND BENEFITS

Babi.Plus BCPAP Kit

Complete kit ready for nCPAP delivery

Includes Bubble PAP Valve and Pressure Relief Manifold

Pressure control settings for easy and secure adjustment of nCPAP

Monitor and adjust fluid level without interrupting BCPAP delivery

Closed system minimizes spillage and atmospheric contamination of H₂O

Integrated oxygen supply tubing with fits all connector for mixed gas delivery

Single patient use, latex free, hypoallergenic material

Babi.Plus nCPAP Cannula Kit

Comfortable method for delivery of non-invasive therapy via nasal prongs for infants weighing ≤ 10 kg

Nasal Prongs made of medical grade silicone to minimize skin irritation

Contoured prong allows better fit without resting on the nasal septum

Single Patient Use kit packaged with four cannulas

Constructed of DEHP-free, latex-free and hypoallergenic material

GiO 1 Pressure Gauge

Real time digital display of pressure measured

Displays peak pressure when Peak and Hold function is active

Auto Power Off preserves battery life

Venti.Plus Pole Clamps

Easily attached to pole with thumb screw mechanism

Available in Single or Dual Clamp

Bubble CPAP (BCPAP) is gaining broad acceptance as it provides a simple method to non-invasively manage the respiratory needs of infants. When combined with early support strategies, BCPAP has proven effective in decreasing the incidence of BPD. Practitioners wanting to use BCPAP as a first choice of therapy find the Babi.Plus system easy to use, safe and effective in performance.

Babi.Plus® BCPAP Kit

Babi.Plus Bubble PAP Valve and Babi.Plus Pressure Relief Manifold



Bubble PAP Valve provides a convenient method of delivering continuous PAP (positive airway pressure) to infants weighing ≤ 10 kg. The patented design delivers accuracy and stability throughout the course of therapy, freeing you to focus on patient care. Adjustable from 1 to 10 cm H₂O, for gas flows from 1 to 12 liters per minute. The Bubble PAP Valve's ergonomic design allows pressure to be easily set without the cumbersome and time consuming task normally associated with bubble devices.

12.5 cm H₂O Pressure Relief Manifold with Gas Supply Tube and Oxygen Sensor Port is for use with patients weighing ≤ 10 kg to prevent breathing circuit pressure from exceeding the preset pressure limit. The inlet adaptor interfaces with all humidifier chambers.



Babi.Plus nCPAP Cannula Kit

Safe, comfortable method to deliver nCPAP and NIV

A single use device used to deliver intermittent or continuous gas flow via a nasal prong to neonates, infants and small pediatric patients weighing ≤ 10 kg.



The contoured nasal prongs are available in 8 sizes to accommodate a wide range of infants; ELBW infant to 10 kg infant. Prongs are made of medical grade silicone to minimize skin irritation and provide stability during non-invasive therapy. The unique cannula body reduces contact with infant's skin while keeping the cannula comfortably in place. An integrated pressure monitoring line allows monitoring of proximal airway pressure during therapy.

GiO™ 1 Pressure Gauge

Simple method for bedside pressure measurement

GiO pressure gauge can be used in a wide range of applications and clinical settings. Battery operated for ease of use and set up. Portable and compact, GiO is designed for both non-invasive and invasive applications at the bedside.




Venti.Plus™ Single, Dual Pole Clamps

Conveniently attach Babi.Plus Bubble PAP Valve

Sturdy pole clamps designed to support the Babi.Plus Bubble PAP Valve and Kit to an equipment pole.

Baby Head Band™ and Circuit Bumpers™

A more comfortable and versatile way to secure nasal prongs

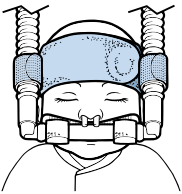

Baby Head Band and Circuit Bumpers secure the nasal interface during oxygen therapy, nCPAP and NIV therapy. The Head Band is made of breathable, non-latex foam to absorb and wick away moisture from the scalp. The soft band conforms to the smallest of heads while maintaining a stable circuit position. Circuit Bumpers secure and stabilize the inspiratory and expiratory circuit limbs to the Band while providing optimal positioning of nasal prongs.

Baby Band™ for Chin and Head


A more comfortable way to close the mouth and position Baby Head Band

Designed to stabilize the position of the Baby Head Band and gently close the mouth during delivery of oxygen therapy, nCPAP or NIV therapy. The soft material of the Band conforms to the chin and head to allow a comfortable fit as needed. The ends of each band allow attachment to the Head Band for best placement.



Baby Nose Bumper™ and RespiraGel Mustache™

For babies who want a pretty nose after therapy

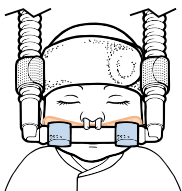

Provides optimal fixation of nasal prongs for infants during nCPAP or NIV therapy. Baby Nose Bumper mustache is made of RespiraGel™, our new "skin friendly" hydrocolloid-based adhesive. Baby Nose Bumper gently holds the nasal interface to the mustache, providing a secure grip and gentle cushion for the nares.

Baby Chin Strap™


For babies who cannot keep their mouths closed

Helps keep small patients' mouths closed during delivery of nasal CPAP and NIV. Offers support with a soft, skin-friendly strap that cups the chin and is secured to the Baby Cap with hook and loop tabs.



Baby Cap and Circuit Bumpers™

For babies who want a soft, comfortable cap

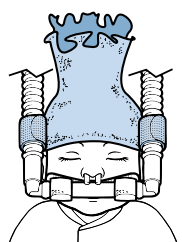

Baby Cap is intended to reduce heat loss and assists in securing the nasal interface during nCPAP or NIV therapy. The soft cotton knit Cap is available in seven sizes and accommodates infants weighing < 10 kgs. Circuit Bumpers gently attach the dual-limb circuit to the Baby Cap, cushioning the infant from the breathing circuit. Baby Cap and Circuit Bumpers allow the circuit and prongs to be maintained in best position for a relaxed fit.

FEATURES AND BENEFITS

Baby Head Band and Circuit Bumpers

Available in four sizes to provide comfortable fit for all infants

Sized by head circumference for easy and quick application

Bumpers designed for secure attachment to Band and minimize pressure points

Open at the crown for fontanel checks and to accommodate scalp vein IVs

Allows caregivers to move the infant without interrupting therapy

Baby Band for Chin and Head

Soft material for patient comfort

Allows caregiver to determine need to close the mouth or position Head Band

Baby Nose Bumper and RespiraGel Mustache

Provides a secure, stable fit of the prongs in the nares

Mustache is constructed of RespiraGel, a latex-free skin-protective hydrocolloid adhesive

Baby Nose Bumper gently grasps the mustache for a secure, stable grip

Baby Chin Strap

Intended for very low birth weight through full-term infants

Assists in keeping the mouth closed to optimize CPAP and NIV therapy

Made of skin-friendly, soft cotton to optimize infant comfort

Single patient use, latex free

Baby Cap and Circuit Bumpers

Fashioned to provide comfortable fit for all infants

Allows caregivers to freely move infant without interrupting CPAP therapy

Crown opening for checking fontanels and accommodating scalp vein IVs

Single patient use, latex free

Ordering Information

The Baby and Babi.Plus products can be used in neonatal critical care units, delivery rooms, pediatric critical care units and special procedure units. Respiralogics products are available from Respiralogics and our specialty respiratory care distributors. For more information visit www.respiralogics.com or email 4info@respiralogics.com.



Babi.Plus Product Line

1002	Babi.Plus Kit BCPAP Kit with Bubble PAP Valve and Pressure Relief Manifold (12/case)
1010	Babi.Plus nCPAP Kit, Small (12/case)
1020	Babi.Plus nCPAP Kit, Large (12/case)
1000	Babi.Plus Bubble PAP Valve (12/case)
1030	Babi.Plus Pressure Relief Manifold (12/case)
B020	Babi.Plus nCPAP Nasal Cannula with Elbow, Size 0 (5/box)
B021	Babi.Plus nCPAP Nasal Cannula with Elbow, Size 1 (5/box)
B022	Babi.Plus nCPAP Nasal Cannula with Elbow, Size 2 (5/box)
B023	Babi.Plus nCPAP Nasal Cannula with Elbow, Size 3 (5/box)
B024	Babi.Plus nCPAP Nasal Cannula with Elbow, Size 4 (5/box)
B025	Babi.Plus nCPAP Nasal Cannula with Elbow, Size 5 (5/box)
B026	Babi.Plus nCPAP Nasal Cannula with Elbow, Size 6 (5/box)
B027	Babi.Plus nCPAP Nasal Cannula with Elbow, Size 7 (5/box)

Venti.Plus Product Line

1005	Venti.Plus Single Pole Clamp (1 each)
1006	Venti.Plus Dual Pole Clamp (1 each)

Baby Product Line

30100	Baby Nose Bumper and RespiraGel Mustache (20/box)
30150	Baby Nose Bumper Replacement (25 bags of 4 each)
30301	Baby Cap and Circuit Bumpers, Size 1 (12/box)
30302	Baby Cap and Circuit Bumpers, Size 2 (12/box)
30303	Baby Cap and Circuit Bumpers, Size 3 (12/box)
30304	Baby Cap and Circuit Bumpers, Size 4 (12/box)
30305	Baby Cap and Circuit Bumpers, Size 5 (12/box)
30306	Baby Cap and Circuit Bumpers, Size 6 (6/box)
30307	Baby Cap and Circuit Bumpers, Size 7 (6/box)
30400	Baby Chin Strap (20/box)
30501	Baby Head Band and Circuit Bumpers, Extra Sm (15/box)
30502	Baby Head Band and Circuit Bumpers, Small (15/box)
30503	Baby Head Band and Circuit Bumpers, Medium (15/box)
30504	Baby Head Band and Circuit Bumpers, Large (15/box)
30510	Baby Band™ for Chin and Head (15 bags of 2 each)
GB10	GiO 1 Pressure Gauge, 0 to 30 cm H ₂ O (1 each)
RL1099	Pressure Line, Transducer Protector Replacement Kit Includes: Hydrophobic Filters (10 each) Pressure Lines (10 each) Straight Connector (1 each)

