Respiragel™ Stoma Cushion

Novel cushion to absorb and redistribute pressure at trachea stoma sites

Stoma sites may become sensitive or compromised due to constant pressure or movement of the tracheostomy tube and flange against tender tissue. The Respiragel Stoma Cushion is intended to cushion the area between the flange and the stoma site, reducing movement and pressure at the site. The contoured surface of the Cushion provides a stable, comfortable interface between the flange and the patient's neck. Respiragel Stoma Cushions are made of medical-grade silicone encased in polyurethane. Hypoallergenic and latex-free, each Cushion is individually packaged and intended for single-patient application. The Respiragel Stoma Cushion is available in two sizes, accommodating children through adults.

The Respiragel Stoma Cushion design allows for placement of the Cushion on the trach tube as early as the time of the tracheostomy procedure. Early use of the Stoma Cushion may assist in reducing irritation and tissue breakdown at the stoma site as well stabilize the tracheostomy tube. Use of the Stoma Cushion may decrease air leaks around the stoma site during trach weaning or during speech therapy by improving the seal between the Cushion and the stoma.

The Respiragel Stoma Cushion may remain in place for up to 30 days and is simply changed at the time of trach change. The Cushion may be used in critical care units, operating rooms, special procedure units as well as alternate care facilities and home care.

The Respiragel Stoma Cushion enhances patient comfort with all styles of tracheostomy tubes. Danny Ties™ may be used in combination with the Respiragel Stoma Cushion to provide a secure fit of the tracheostomy tube.


Ordering Information

50005	Respiragel Stoma Cushion, 5 mm	10/box
50005C	Respiragel Stoma Cushion, 5 mm	40/case
50010	Respiragel Stoma Cushion, 10 mm	10/box
50010C	Respiragel Stoma Cushion, 10 mm	40/case

Companion Products

DWL105001	Danny Ties, Small
DWL105002	Danny Ties, Medium
DWL105003	Danny Ties, Large
DWL105005	Jimmy Proof Strap™

Respiralogics products are available from Respiralogics and our specialty respiratory care distributors. For our latest product information visit www.respiralogics.com or email us at 4info@respiralogics.com.

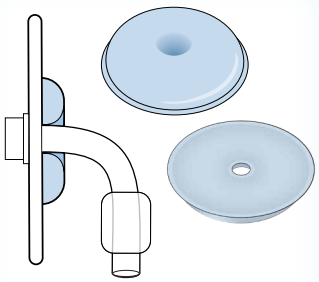

FEATURES AND BENEFITS

Individually packaged, available in two sizes and 5.0 mm and 10.0 mm thicknesses, accommodating children through adults

The Stoma Cushion seals around the tracheostomy tube up to 16.0 mm outside diameter

Contoured surface rests against the neck and reduces flange contact with stoma site

Hypo-allergenic, latex free, medical-grade silicone Cushion encased in polyurethane

Designed to be placed on the tracheostomy tube prior to placement of the trach tube and changed at time of trach change

30-day use and easily cleaned

Respiralogics is a Global Respiratory Solutions, Inc. Company
3545 Airway Drive, Suite 104 • Reno, NV 89511 • 775 954 0160 • www.respiralogics.com

© 2026 Respiralogics. All rights reserved. Respiralogics, the Respiralogics logo, and RespiraGel are trademarks of Global Respiratory Solutions, Inc. Danny Ties and Jimmy Proof Straps are trademarks of Leckie Medical Products, Inc. C50005-1