Danny Ties™

Unique tube holder with quilted design minimizes irritation and stabilizes trach

Danny Ties are unique tracheostomy tube holders with a softer and more comfortable fit around the neck for patients of all ages. The patient with a tracheostomy needs to have a tube holder that securely holds the artificial airway in place to prevent accidental decannulation. As important, the tube holder needs to provide a soft, comfortable fit about the neck while minimizing skin irritation under the collar. Danny Ties are made of soft, absorbent cotton that lays smooth at the edges of the collar, minimizes skin irritation and reduces skin breakdown under the collar. The patent-pending design of the Danny Ties evenly distributes the quilted collar around the neck to minimize pressure points on the skin.


Danny Ties are easy to apply with tapered ends on the collar straps. The ends thread easily through the eyelets of the tracheostomy tube, allowing for quick application and changes of the collar position on the smallest of infants and the largest adult patient. The collar is engineered to hold its shape—not fold in half around the neck when the patient turns—and does not stretch when it absorbs moisture. If soiled, Danny Ties may be washed, dried and reused for the same patient per institutional policy.

The Respiragel™ Stoma Cushion may be used in combination with the Danny Ties to provide a secure fit of tracheostomy tube. Jimmy Proof Straps™ may be used with the Danny Ties to reduce risk of accidental disconnects from the ventilator circuit.

Danny Ties are developed with special care by a dad wanting to "make a difference" for his son and change the quality of life for all tracheostomy patients. His engineering expertise and commitment to finding a better solution is found in the Danny Ties.

Sizing and Fitting Guidelines

Product Code	Size	Fits Neck
DWL105001	Small	6–10 inch
DWL105002	Medium	9–20 inch
DWL105003	Large	19–25 inch

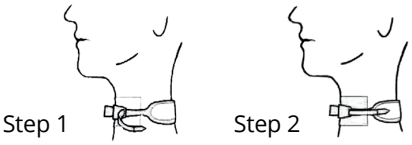

Ordering Information

DWL105001	Danny Ties, Small	100/case
DWL105001B	Danny Ties, Small	10/box
DWL105002	Danny Ties, Medium	100/case
DWL105002B	Danny Ties, Medium	10/box
DWL105003	Danny Ties, Large	100/case
DWL105003B	Danny Ties, Large	10/box

Companion Products

DWL105005	Jimmy Proof Straps™
50005	Respiragel Stoma Cushion, 5 mm
50010	Respiragel Stoma Cushion, 10 mm



FEATURES AND BENEFITS

Danny Ties stabilize the tracheostomy tube against the trachea with one-piece design

Quilted collar with reinforced stitching provides a softer, more comfortable fit about the neck

Unique tapered straps easily thread the ties through the tracheostomy flange

Absorbent cotton collar wicks moisture away. Helps to minimize irritation and reduce skin breakdown around the neck

Reinforced stitching at the ends of the tie strap to ensure ease of threading and repositioning of the collar

Available in three sizes to accommodate small infants through adults

Single patient use

Danny Ties may be gently washed, dried and reused for the same patient

Respiralogics products are available from Respiralogics and our specialty respiratory care distributors. For our latest product information visit www.respiralogics.com or email us at 4info@respiralogics.com.

Respiralogics is a Global Respiratory Solutions, Inc. Company
3545 Airway Drive, Suite 104 • Reno, NV 89511 • 775 954 0160 • www.respiralogics.com

© 2026 Respiralogics. All rights reserved. Respiralogics, the Respiralogics logo, and Respiragel are trademarks of Global Respiratory Solutions, Inc. Danny Ties and Jimmy Proof Straps are trademarks of Leckie Medical Products, Inc. CDWL105001.6